

LM-6800

Product Massage Chair

Model LM6800

Distributed by AJX Distribution LLC

7001 Village Dr #250

Buena Park, CA 90621


Tel 1-855-646-9500

Email INFO@KAHUNACHAIR.COM

KAHUNA MASSAGE CHAIR®

LM-6800

Massage Chair Manual

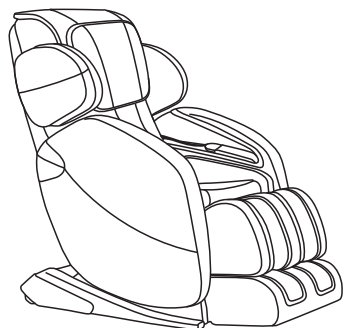

INSTRUCTION MANUAL
MODEL:LM6800


Please read this user manual thoroughly prior to operating the massage chair. Store this booklet for future reference.



KAHUNA  CHAIR


MAINTENANCE

CONTENTS

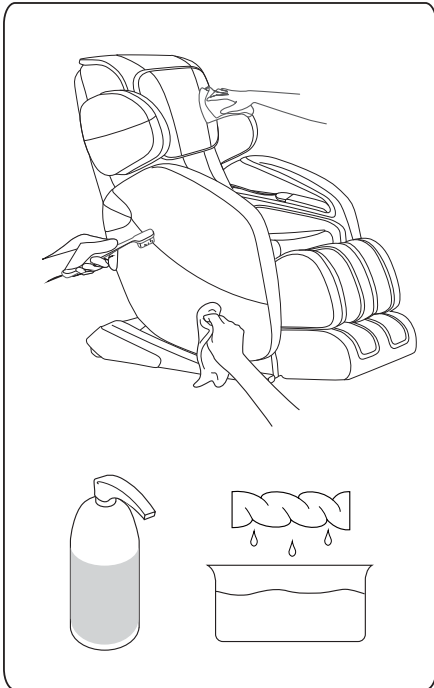

Use dry cloth to clean

Turn off and pull out the plug before cleaning the massage chair, do not touch the plug with wet hands.

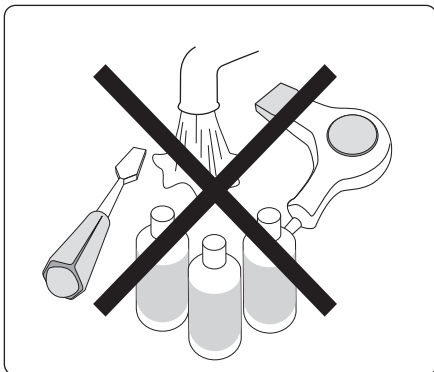

Prohibit using liquid cleaning products like water, alcohol, or other cleaning product to maintain a good quality.

Please read and follow all safety and operating instructions listed in this user manual prior to use. Basic precautions should always be used when operating an electrical appliance.

Safety Instructions	2
Product Specifications	4
Product Components	5
Chair Installation	7
Remote Controller	13
Operation Guide	14
Troubleshooting	19
Maintenance	20

SAFETY INSTRUCTION

Warning

The power plug must be inserted properly to the outlet. If the connection of the equipment-grounding conductor is not connected properly it can result in the risk of electric shock. Check with a qualified electrician or service technician, if you are in doubt as to whether the product is properly grounded. Do not modify the plug that was provided with the product, if the outlet is not suitable. Please make sure that the outlet is properly installed by a qualified electrician. When using the massage chair, basic precautions should always be followed to reduce the risk of electric shock, burns, fire or injury.

1. Please make sure that the power plug of massage chair is properly inserted in a proper electrical outlet.
2. Make sure the power plug is properly inserted into the proper electrical outlet to reduce the risk of short-circuit and fire.
3. Always unplug the massage chair from the electrical outlet immediately after use, before cleaning, etc.
4. Do not use the massage chair in places where the humidity is high.
5. Do not use the massage chair outdoors and avoid place with direct sunlight.
6. If the power turns off during use make sure to switch the power OFF position and pull the power plug out from the socket to prevent damage from occurring when the power is restored.
7. Do not operate the massage chair if the power cord is wet or damaged.
8. Do not plug/unplug the power cord from the socket with wet hands.
9. Supervision is necessary when the massage chair is being used by on or near the children or disabled persons.
10. Use the massage chair only for its intended use as described in this user manual.
11. Do not stand, sit or jump in the calf rest, back rest and armrest when the massage chair is in use to prevent any injury.
12. Do not fall asleep during the massage session.
13. Always keep the massage chair clean and do not insert any object into the gaps and around the massage chair.
14. Do not move the massage chair by pulling the power supply cord.
15. Make sure that there is no person (children in particular) pet or obstacles underneath the massage chair prior to lowering the leg rest/back rest.

TROUBLE SHOOTING


If the sounds or movements are not normal, please turn off the massage chair and remove the plug from the electric outlet. Please contact the customer service department for further assistance.

855.646.9500 / 1.714.388.2625, service@kahunachair.com

Problem	Solution
The massage chair is connected to power but is not functioning.	Check to see if the power plug is inserted or if the power outlet is functioning properly.
The remote controller is responding, but the massage chair is not functioning.	The massage chair may have automatically shut down due to overheating protection devices. Please turn off the power and allow the chair to cool off for 30 to 50 minutes.
Both the remote controller and massage chair are not responding.	The massage chair automatically turns into "full stop" mode when a massage program is completed. Please press the power button again to turn of the chair.


OPERATION GUIDE

3 After use

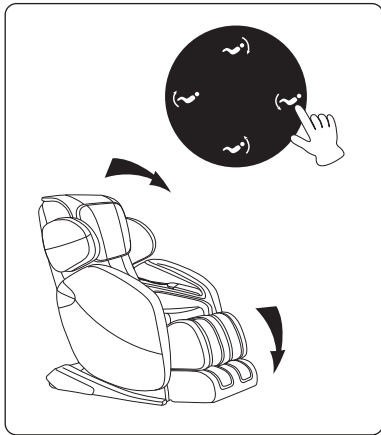


1. After the massage program is completed return the massage chair to upright position by turning off the massage chair.

Please check the surrounding prior to adjust the recline angle, paying special attention to children or pets that may be in close proximity.

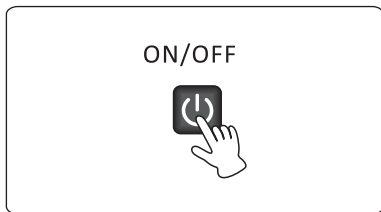

2. Press the power button to power off the massage chair and return remote controller to the storage area on arm pocket.

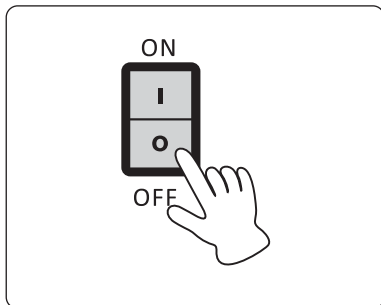

3. Turn off main power switch when you are not using the chair for a long period of time. Remove Plug from the power outlet to avoid accidents by kids and pets.

4. Do not operate with wet hands

TROUBLE SHOOTING

16. Always check the surface of the massage chair if it is torn or ripped before usage, if you find any unfamiliar rips in any area of the cloth in the back rest, do not use the massage chair. Kindly contact the service center for repair.
17. This appliance is designed for household use only.
18. For those suffering from illness, please consult your doctor for advice prior to use.
19. The following conditions should consult doctors before use:
 - (a) Suffering malignant tumor
 - (b) People with heart disease
 - (c) Women in pregnancy, menstruation and puerperal fever
 - (d) People who are medicated or under medication
 - (e) People with osteoporosis
 - (f) People with back bended exceptionally by the accident or disease
20. If you experience pain during usage, stop using the massage chair immediately. If any injuries occur, seek medical treatment immediately.
21. Do not massage the swollen or inflamed parts of the body.
22. This product is non-professional appliance designed to provide a comfortable massage. Do not use this product as substitute for medical treatment.
23. Each massage session should not exceed 30 minutes. Massage the targeted area only for 15 minutes, over massaging can lead to over heating, thus giving the appliance a shorter lifespan.
24. Do not operate where aerosol (spray) products are used or where oxygen is being administered.
25. The massage chair is not suitable for children age under 12 and below.
26. Please make sure that the massage chair is placed on a flat level ground.

NOTE

This appliance is not intended to use by people (including children) to reduce physical sensory or mental capabilities or lack of experience and knowledge, unless they have been given supervision or instruction concerning to use this appliance by a responsible person for their safety.


Children must be supervised to make sure that they do not play the appliance. If the supply cord is damaged it must be replaced by a special cord or assembled available from the manufacturer or its service agent. This appliance heats up during operation, people sensitive to heat must take extra caution when using this appliance.

PRODUCT SPECIFICATIONS

Description	Specifications
Product Name	Kahuna Massage Chair
Model No.	LM 6800
Rated Voltage	110-120V~
Rated Frequency	50/60Hz
Rated Power Consumption	230 W Max.
Dimension (LxWxH) Inches	46 x 31 x 48 inches
Weight (LBS)	Net Weight: 200 lbs Gross Weight :220 lbs

OPERATION GUIDE

2 Display Diagram

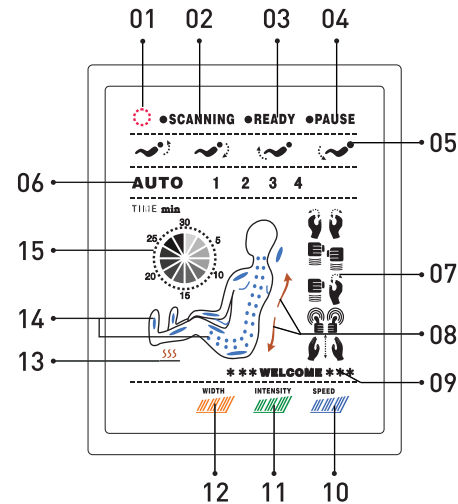

1. Massage chair is powered on.
2. Massage chair is scanning the back and preparing for initiation.
3. Scanning is complete and ready for massage programs.
4. Pause button has been pressed and chair is temporarily stopped.
5. Reclining position adjustment buttons are pressed and the back-rest or leg rest is moving.
6. Display which automatic program is activated.
7. Display which massage technique is activated. Top to Bottom: Kneading, Tapping, Dual, Rolling, Shiatsu.


8. Display the direction of the massage rollers, moving upwards or downwards.
9. Greetings message when preparing to start.
10. Shows the massage speed levels, from up to 3 levels.
11. Shows the massage speed level, from one bar to three bars illuminated.
12. Shows massage roller width, from narrow to medium, to wide.
13. Heating therapy function is on
14. Display the location of massage rollers and activate air bags.
15. Display current time remaining in current massage program.

OPERATION GUIDE

Annular functional area



Air pressure massage area button. Each time you press this button, namely to lit the air pressure massage area on the screen. Press the button continued, there will be different parts of the combination of massage.

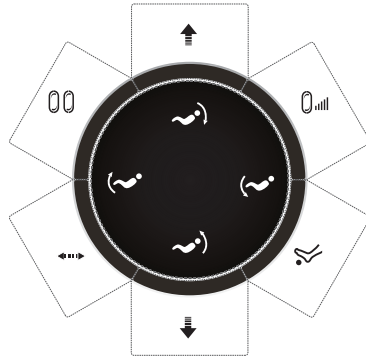

This button allows for you to control the intensity of the air bag on 3 different intensity levels.



Press this button to turn on/off the foot roller function.



Roller location adjustment: Hold the up or down arrow to move the massage rollers to the desired location.



Width: Press to adjust massage roller widths from narrow to medium to wide. Width adjustment is not available in kneading mode.



Press the button and hold on to adjust the foot rest independently.




Press the button and hold on to adjust the massaging position of back-rest and leg rest.

P.S: The position adjustment button will not work when the chair in the yoga stretch program.


PRODUCT COMPONENTS

Main Body Diagram (Front)

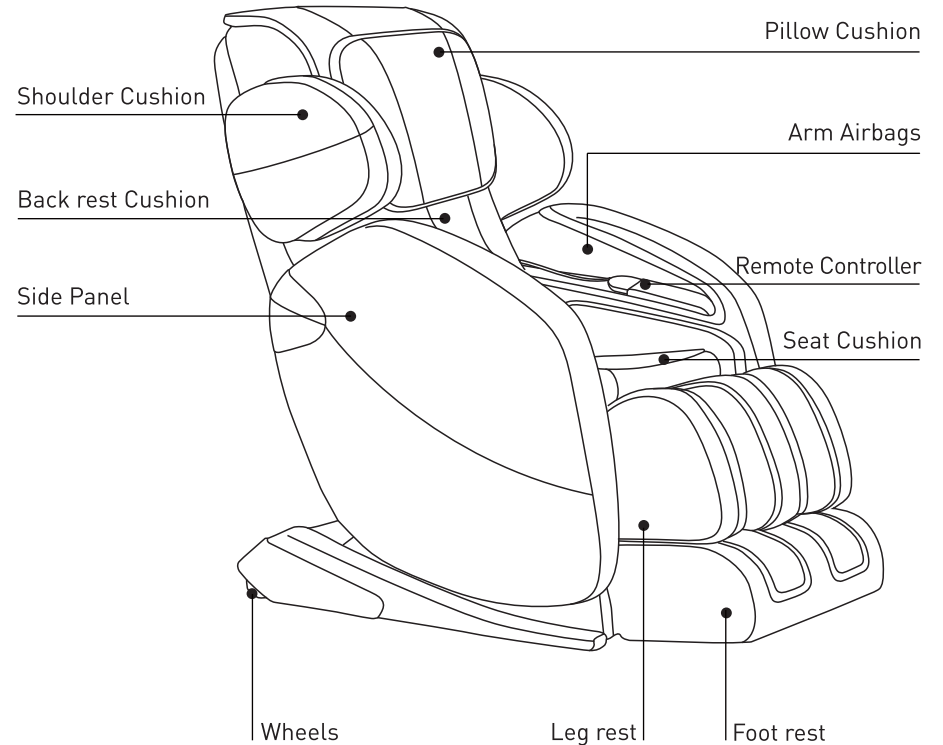

PRODUCT COMPONENTS

Main Body Diagram (Rear)

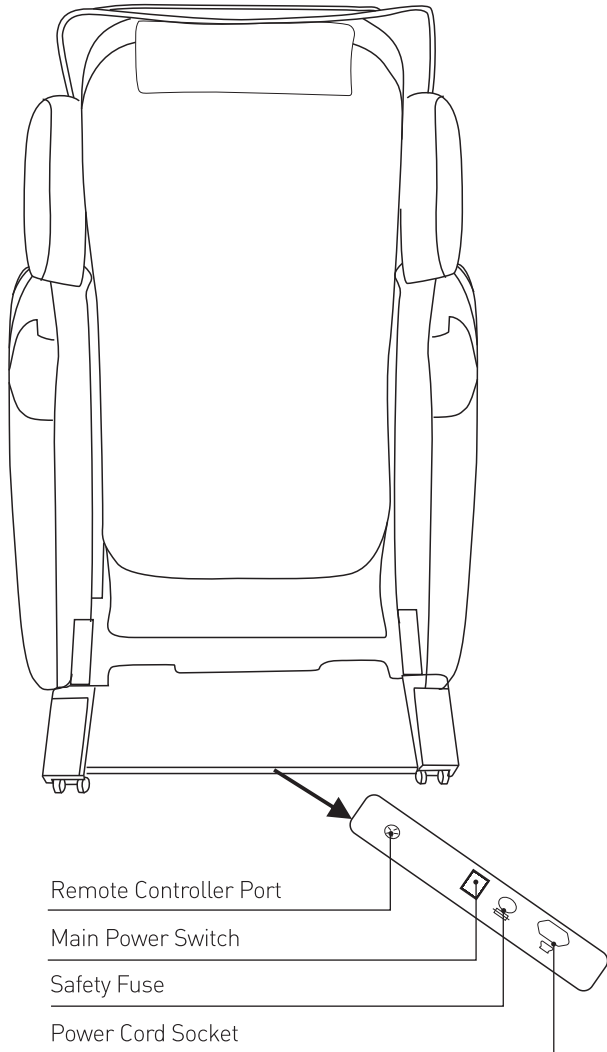

OPERATION GUIDE

This chair design with 2 sets of auto massage program



Yoga stretching mode stretches your full body, LM-6800 Signature auto mode.



SV-Chiro Mode helps your back like chiropractic mode. It is designed to release subluxations in the spine to make sure the body and the nervous system are operating at their best possible level.

Main Body Diagram (Rear)



Manual massage technique button,
 Kneading → Tapping → Dual → Rolling → Shiatsu → Off



Speed for Tapping, Kneading, Dual action and Shiatsu speed adjustment up to 3 levels.



“Area adjustment button” In manual mode, press to select between specific areas in the following order:

Back → Waist → Seat → Back & Waist → Waist & Seat → The whole body



“Fixed Spot Massage” Press to stop the rollers at a particular area for a fixed point focal massage in manual mode.



Heating Therapy Switch on/off for back heating function.

OPERATION GUIDE

1. Function description of Controller buttons

Start or Stop massage chair

Press this button, massage chair will be ON, screen will be ON, User can choose auto or manual function to massage.

Press this button again, massage chair will be off, massage stopped, massage chair will come back to the original position, all buttons are inactive. When the chair returns to the original position the chair will completely shut off.

Press to temporarily pause the massage, chair will end idle mode. Press again to resume massage.

Press the button once, massage back and leg position to zero G1, press twice, the position to zero G2 twice, press third times, the position to zero G3, press fourth times, the massage chair return to original position.

Zero G1 → Zero G2 → Zero G3 → Original Position

PS: The position adjustment button have no function when the chair in the stretch program.


Auto mode

Intelligent Automatic Programs: Press once to initiate the first of four intelligent automatic programs. Press repeatedly to switch between the programs, which are sequenced as follow:

Deep Tissue → Stretching → Relax → Rejuvenate → Off

CHAIR INSTALLATION

1. Remove the package and take out the spare parts

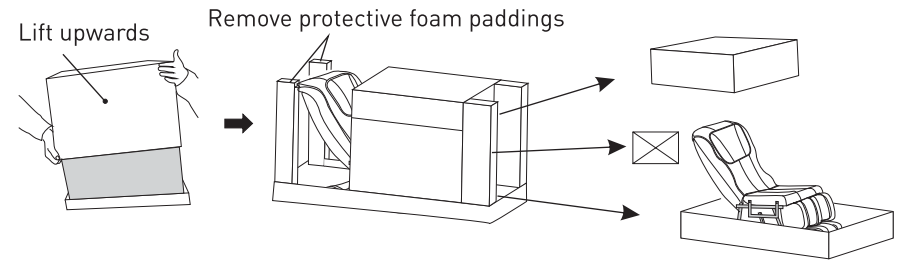


Following contents from

STEP1 Remove the cartons.

STEP2 Install the pillow, back rest cushion, cushion and spare parts.

STEP3 The spare parts include: remote control, user manual and power cord.

2. Install the side board

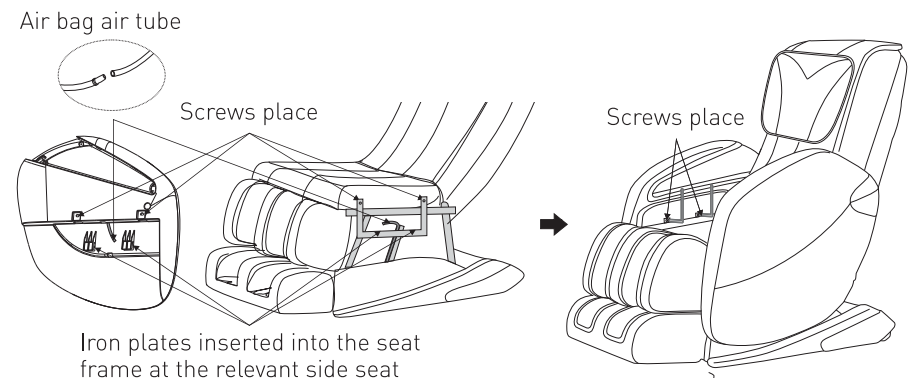


STEP1 Connecting the side plate main body frame, must insert the end in order to prevent air leakage.

STEP2 Aligning the side plate hole and the massage chair hole, then inserted the plate into the corresponding side mount, use the screw driver lock-in screw fixed. The other side completed according to the above mentioned.

STEP3 Please check the air hose is connected properly to avoid air leakage.


CHAIR INSTALLATION

3. Install the shoulder airbag

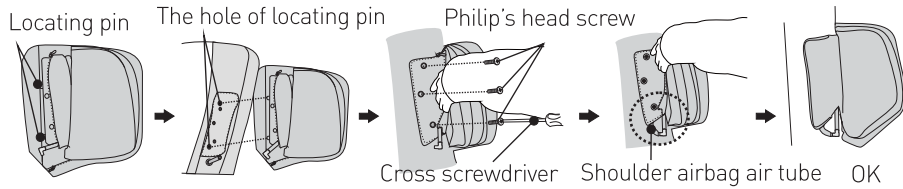


- STEP1** Find the 2 locating pin and screw holes
- STEP2** Using the locating pin insert into the chair back rest, then a cross screw driver to fix the Philip's head screw.
- STEP3** Connect the tube with tube connector on the chair.
- STEP4** Zip up the opening after check each parts are secured.
- STEP5** Another side will follow the same steps.
- STEP6** Double check the installation steps, smooth out any creases in the cushions.

4. Place the massage chair in an open area

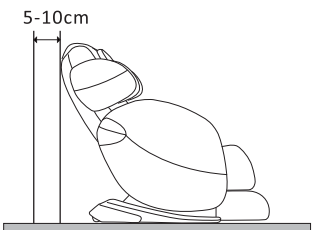

- Allow sufficient space around the chair for recline purposes.
- 5-10cm distance for the back side and 50cm in front of the foot rest.

5. Moving the massage chair

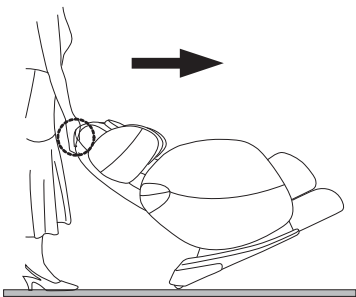

- As illustrated, push the back rest downwards and tilt chair on to its wheels at 45° angle. Slowly push and maneuver the massage chair to desired location. Lower massage chair slowly and smoothly until it is fully flat on the ground.
- CAUTION:**
- *Do not move massage chair when occupied or in use.
 - *Moving the massage chair may cause damage to the floors. Floor protection is recommended (e.g. rugs or carpet).
 - *Do not drop or suddenly let go of the chair when moving as it may damage the chair.
 - *Two persons are recommended when moving and extra care when handling the massage chair.

REMOTE CONTROLLER

1. Remote Controller



- 1. On/off Massage Chair
- 2. Zero Gravity Button
- 3. Auto Mode
- 4. Yoga Button
- 5. Mode
- 6. Heating Button
- 7. Location Adjust Button
- 8. Air Massage Intensity
- 9. Recline Adjust
- 10. Foot roller Button
- 11. Display
- 12. Modes of Manual Program
- 13. Speed Adjustment
- 14. Area Adjustment Button.
- 15. Fixed Spot Massage
- 16. Pause/emergency Stop
- 17. Air Pressure Button
- 18. Roller Width

CHAIR INSTALLATION

Tube Connection Check List

2. Correct tube location connection



Please change the tube connection like the pictures shown above. If tube connection position is correct you don't need to change it.

CHAIR INSTALLATION

Install Guide

Components



BODY



SLIDE PANEL X 2



SLIDE PANEL LEFT X 1 / RIGHT X 1




SHOULDER PANEL
SCREWS X 3 EACH



SIDE PANEL SCREWS
X 2 EACH

Both side panels installation x 2

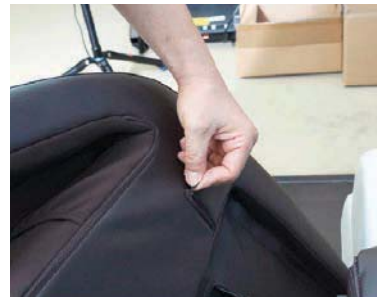

Open the Zip on side panel



Connect Airtubes



Close the zipper



Tighten the screw make sure holes
are stright to screws

CHAIR INSTALLATION

Install Guide

Shoulder panel installation x 2



Open the Zip on side panel



Connect Airtubes



Tighten the screws
Make sure holes are
straight to screws



Connect remote controller



Arrange seat cushion neatly



Close the zipper



Connect remote controller
cable back of the chair




Connect main power cable

CHAIR INSTALLATION

Tube Connection Check List

1. Incorrect tube location connection

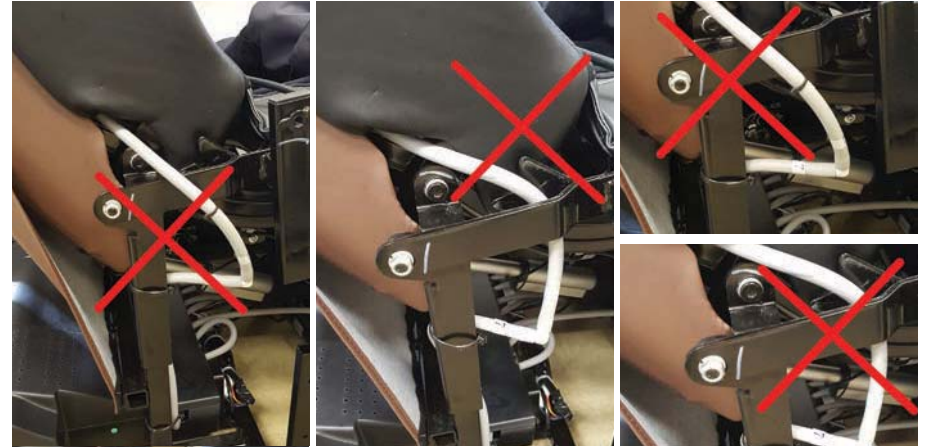

Tube must come out behind the visible frame, please see the image below